



1 Create Your Account

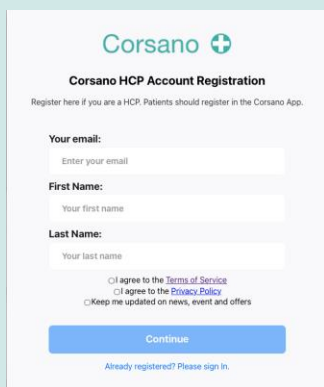


1) Go to <https://portal.corsano.com/#/registration-hcp-admin> and enter your details. The email address will become your username.

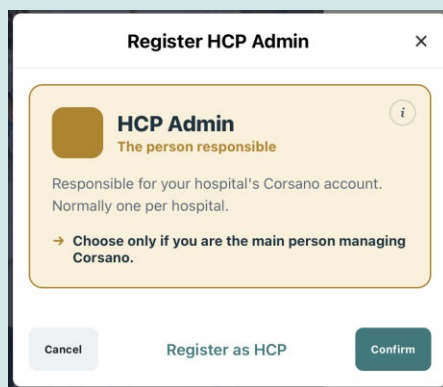


2) Choose HCP Admin as a role and click on confirm

Getting Started – Portal Creation

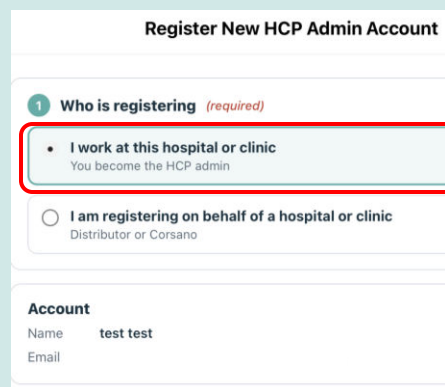


2) Choose "I work at this hospital or clinic" and continue the registration.

2 HCP Admin Account Information

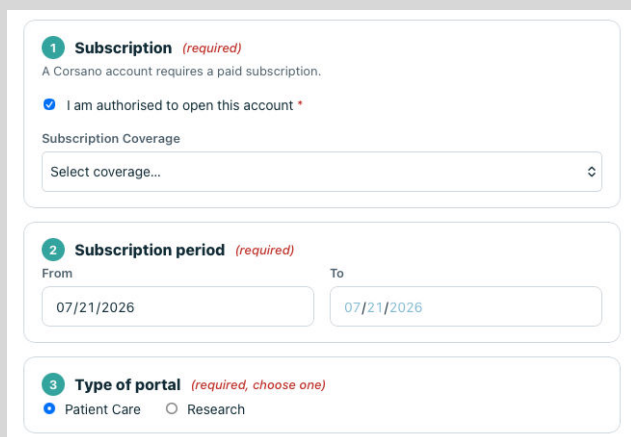


1) Enter the requested information according to your agreement with Corsano or its distributor.

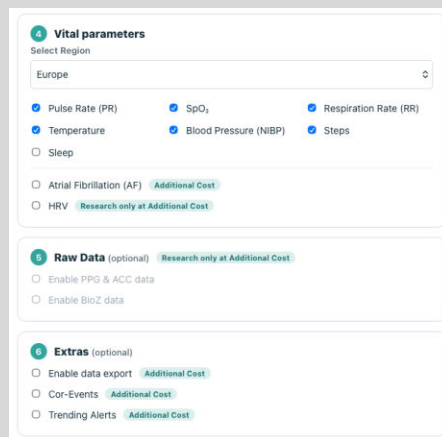


2) Select the vital signs and additional features included in your agreement.

3 Complete the Registration

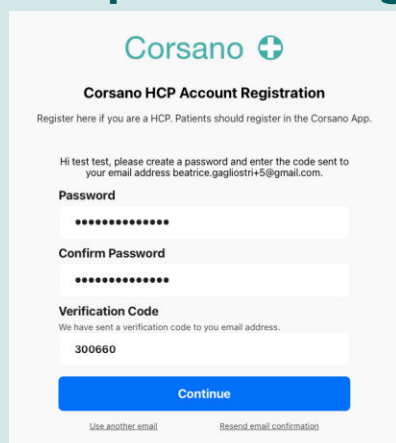

1) Enter the code sent to your email and create a password.

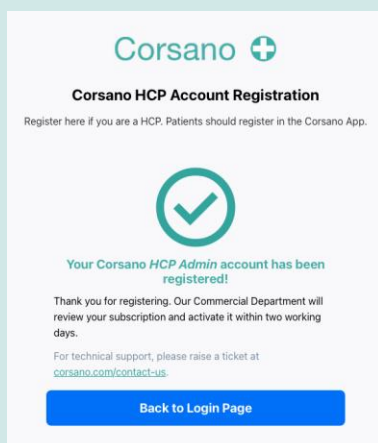

2) Registration complete! You'll receive an email once your account is approved.

Corsano +

Tutorials

www.corsano.com/patient

Support

www.corsano.com/support


**Follow Instructions
in the Corsano App**

Corsano Health B.V.


Waldorpstraat 5
2521 CA Den Haag
The Netherlands

Getting Started – Login & Voucher

1 Login in the Portal

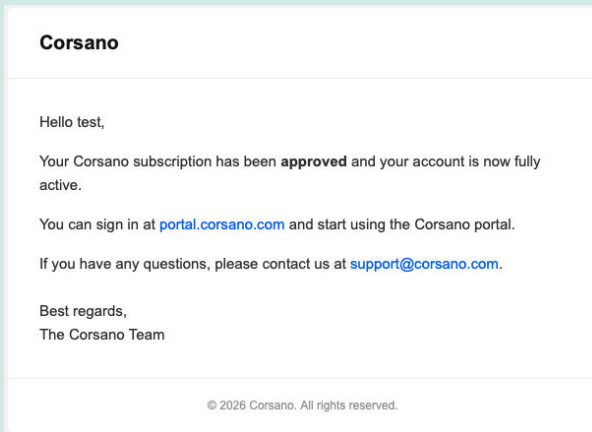

1) You will receive an email once the authorization is granted.

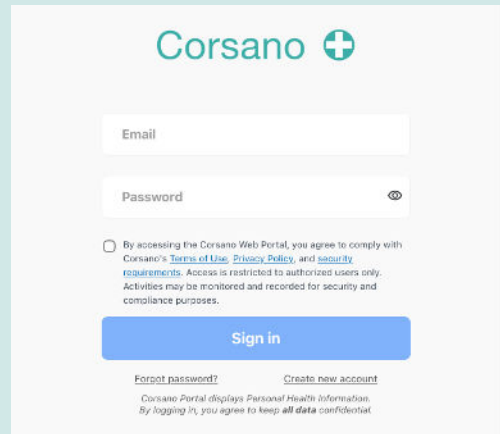


2) Log in in the portal using your registered email and password.

2 Create a Study/Department and a Voucher

Departments


Add Department

1) Create a Study/Department: Click Add Department, enter the study/department name, and click Add.

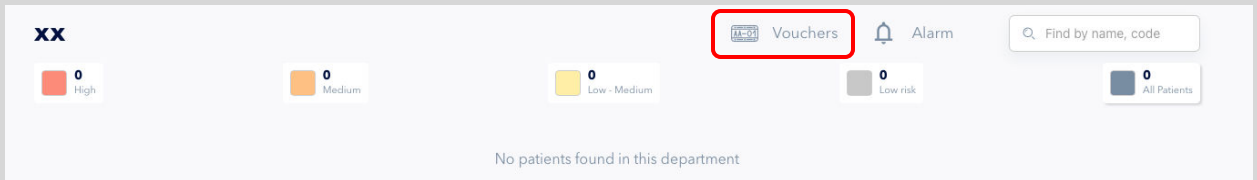


2) Create a Patient (i.e. Voucher) in the study. Open the study/department, go to Vouchers, select the number of vouchers, and click Add Users.

3 Print the Voucher QR Code

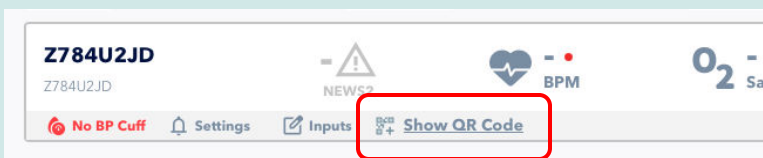

1) Patient created! To view the QR code, click the QR code icon next to the patient.

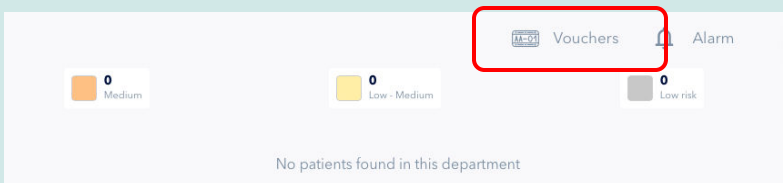

2) You can also click on the Vouchers button and select Patient Letter to print it. Then follow the Getting Started – Pair the device guide to pair the device.

4 Connect

To receive updates regarding all the Corsano Systems, register in Corsano Connect:
<https://connect.corsano.com/register>

Corsano +

Tutorials

www.corsano.com/patient

Support

www.corsano.com/support

Corsano Health B.V.

Waldorpstraat 5
2521 CA Den Haag
The Netherlands

Getting Started – Pair the bracelet


1 Register



1) Open the Camera App on your mobile




2) Take your letter and scan the QR-Code

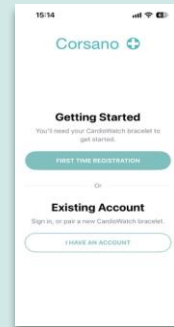

3) Install the App Then select [First Time Registration]




4) Scan the QR Code again with the App

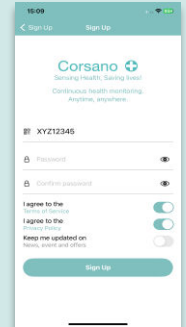

5) Set your own Password and store in safe place!

2 Pair Bracelet & BP Cuff

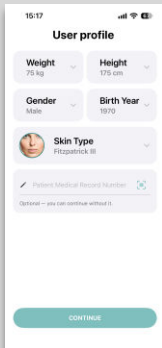

1) Fill in your user profile information



2) Scan the QR code on the case of your bracelet to pair your bracelet



3) Press and hold the Start button of the BP Cuff Monitor for 5 sec to pair your BP Cuff



4) Wear the bracelet on the right arm and the BP Cuff on the left arm



5) Press the Start button of the BP Cuff Monitor to start the BP Calibration

3 Wearing Tension & Position



1) Tighten the strap so that it is secure against the skin, the green sensor light should not be visible. The Bracelet should be snug but comfortable



2) Position: Wear the Bracelet approx. 1 cm below the wrist bone

Corsano +

Tutorials

www.corsano.com/patient

Support

www.corsano.com/support

Follow Instructions in the Corsano App

Corsano Health B.V.

Waldorpstraat 5
2521 CA Den Haag
The Netherlands